

# F&P AIRVO™ 2

Delivering Nasal High Flow with Optiflow™




## F&P AIRVO 2

humidifier with integrated flow generator,  
delivering Nasal High Flow with Optiflow™



# Introducing the F&P AIRVO 2, delivering respiratory support and airway hydration with a focus on comfort and ease of use.

The F&P AIRVO 2 has been designed specifically to be an integrated solution for delivering Nasal High Flow with Optiflow™. Fisher & Paykel Healthcare's leading humidification technology allows the AIRVO 2 to comfortably provide high flows of air/oxygen mixtures to spontaneously breathing patients, through the unique Optiflow™ nasal interface.

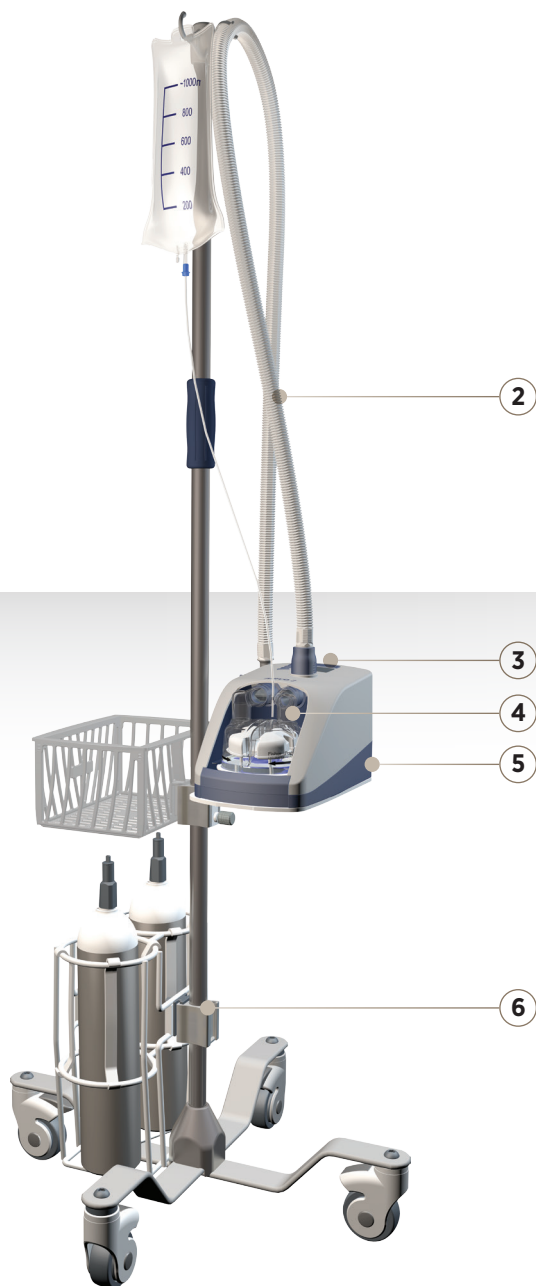

### 1. Optiflow™ Nasal Interface

- Soft, flexible prongs
- Wide-bore design reduces gas jetting
- Prongs contoured to the patient's nose

### 2. Heated breathing tube

- Dual spiral heater wires and unique integrated temperature sensor
- No separate temperature probes or heater-wire adapters required

### 3. Designed for simple setup, use and cleaning

- Helpful onscreen animations assist with setup and troubleshooting

### 4. Adjustable temperature and flow settings

- Three temperature settings (37, 34, 31 °C) help achieve comfort and compliance
- Integrated flow generator delivers a wide flow range (2 - 60 L/min) for both pediatric and adult patients - no wall air supply required

### 5. Supplementary oxygen (when required)

- Oxygen can be added from a wall supply or cylinder
- Inbuilt ultrasonic oxygen analyzer requires no calibration, service or replacement

### 6. An integrated system for delivering Nasal High Flow with Optiflow™

- The F&P AIRVO 2, with its wide range of accessories, provides versatility, mobility and convenience.